



Sterile Container System & Baskets



Classic & Bio-Barrier Model

We know the value of life!



Thempson

Surgical Instruments

Surgical innovation Sytems

From Access to Repair

Address:

Veko Giz Plaza, Meydan Sokak No.3/43, Kat.13 Maslak, Istanbul,Turkey

Phone: +90 212 705 36 16, Fax: +90 212 705 36 36

info@thempson.com

www.thempson.com


ISO 13485:2003

ISO 9001 : 2008



- According to the production programme Thempson`s rights is reserved in production to improve, modify or delete the the products declared in this catalog.

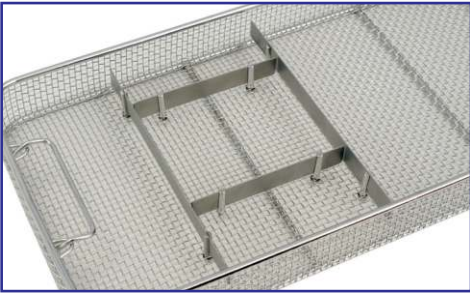- The weight information given in our catalog is approximate value.



The holding pins and plates for separation wire baskets or desire



The silicon mat, for placing the micro instruments.



Silicone handle has available in six different colors and detachable by user.



The drape retainer.



Spiral holder



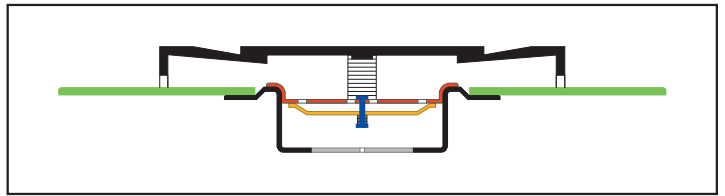
Wire bottom, balanced placement and steam circulation safety.




The fixation clamps for the cylindirical instruments in different sizes.




The method of mounting - de-mounting the safety lid of the fix model container.



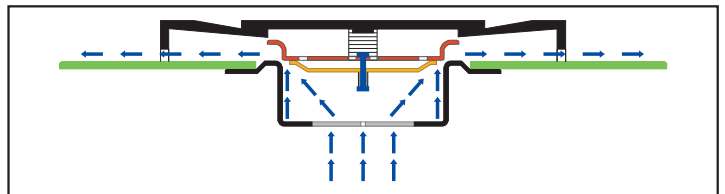
1. Before sterilization the valves are in relax position.



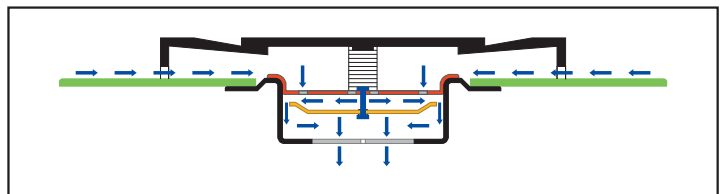
2. During the sterilization by effect of vacuum the valves are opened completely and the dry air goes out of the container.



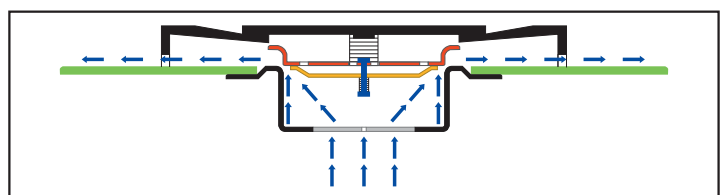3. At the end the sterilization cycle the valves are opened in the reverse way by effect of vacuum and the sterilization cycle is completed.



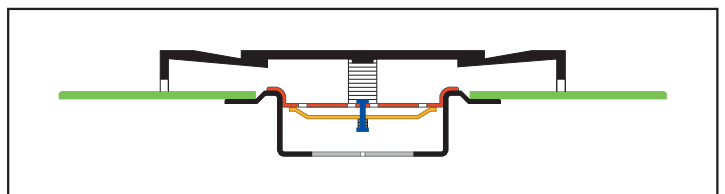4. At the end of the sterilization cycle the sterilizer equates the pressure level and the valves of the Barrier Model Container are closed.



5. After sterilization the valves are in relax position.

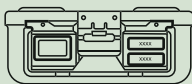
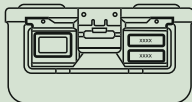
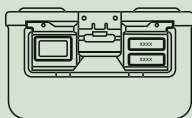
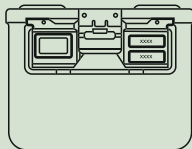
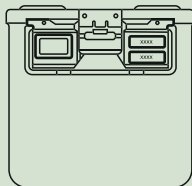



Sterile Container System




1/1 Full Size

Classic Model

Specifications Sizes	Lid (Perforated)	Lid (Perforated)	Safety Lid	Safety Lid
	Box (Non-Perforated)	Box (Perforated)	Box (Non-Perforated)	Box (Perforated)
580 X 280 X 100 mm 	10-010  10-011  10-012  10-013  10-014  10-015 	10-040  10-041  10-042  10-043  10-044  10-045 	10-150  10-151  10-152  10-153  10-154  10-155 	10-180  10-181  10-182  10-183  10-184  10-185 
580 X 280 X 135 mm 	10-016  10-017  10-018  10-019  10-020  10-021 	10-046  10-047  10-048  10-049  10-050  10-051 	10-156  10-157  10-158  10-159  10-160  10-161 	10-186  10-187  10-188  10-189  10-190  10-191 
580 X 280 X 150 mm 	10-022  10-023  10-024  10-025  10-026  10-027 	10-052  10-053  10-054  10-055  10-056  10-057 	10-162  10-163  10-164  10-165  10-166  10-167 	10-192  10-193  10-194  10-195  10-196  10-197 
580 X 280 X 200 mm 	10-028  10-029  10-030  10-031  10-032  10-033 	10-058  10-059  10-060  10-061  10-062  10-063 	10-168  10-169  10-170  10-171  10-172  10-173 	10-198  10-199  10-200  10-201  10-202  10-203 
580 X 280 X 260 mm 	10-034  10-035  10-036  10-037  10-038  10-039 	10-064  10-065  10-066  10-067  10-068  10-069 	10-174  10-175  10-176  10-177  10-178  10-179 	10-204  10-205  10-206  10-207  10-208  10-209 
Ps: All container measurements cover the outer size of the lids.				
				

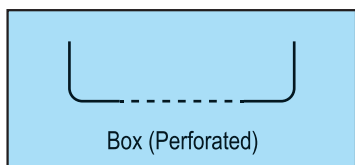

CODE	OUTSIDE	INSIDE	WEIGHT
10-100	595 X 275 X 95 mm	543 X 260 X 75 mm	2.15 kg
10-101	595 X 275 X 116 mm	543 X 260 X 96 mm	2.30 kg
10-102	595 X 275 X 135 mm	543 X 260 X 115 mm	2.40 kg
10-103	595 X 275 X 185 mm	543 X 260 X 165 mm	2.75 kg
10-104	595 X 275 X 243 mm	543 X 260 X 223 mm	3.25 kg

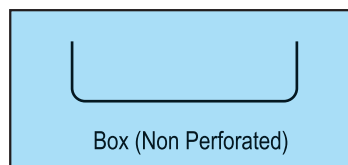

CODE	OUTSIDE	INSIDE	WEIGHT
10-105	595 X 275 X 95 mm	543 X 260 X 90 mm	1.70 kg
10-106	595 X 275 X 116 mm	543 X 260 X 111 mm	2.05 kg
10-107	595 X 275 X 135 mm	543 X 260 X 130 mm	2.15 kg
10-108	595 X 275 X 185 mm	543 X 260 X 180 mm	2.60 kg
10-109	595 X 275 X 243 mm	543 X 260 X 238 mm	3.10 kg



Perforated Lids



Safety Lids

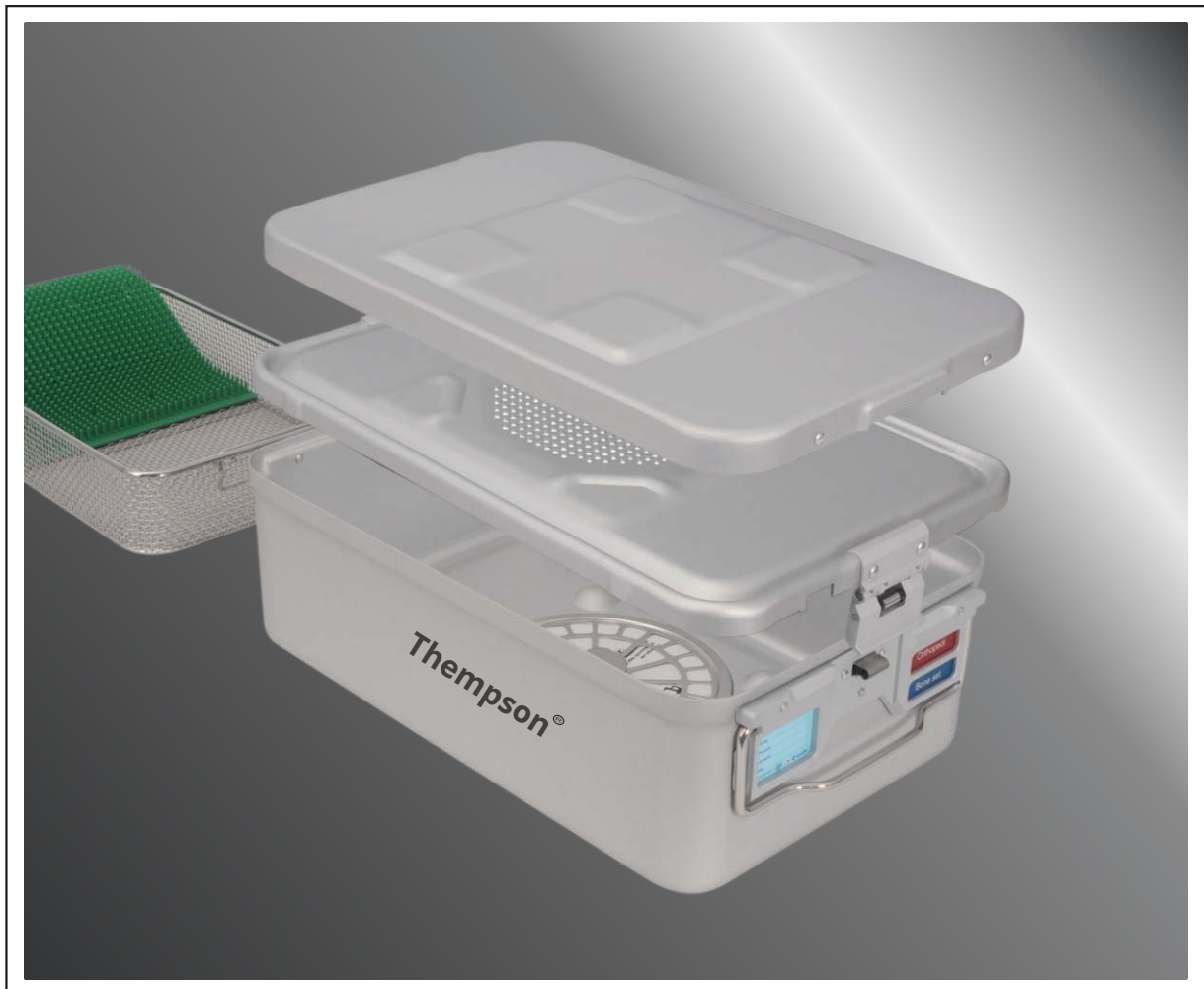


Removal Lids (Non Perforated)









Sterile Container System




3/4 Quarter Size

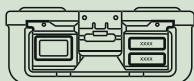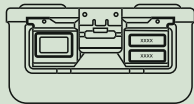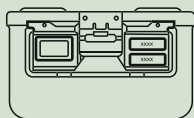Classic Model

Specifications Sizes	Lid (Perforated)	Lid (Perforated)	Safety Lid	Safety Lid
	Box (Non-Perforated)	Box (Perforated)	Box (Non-Perforated)	Box (Perforated)
465 X 280 X 100 mm 	10-240  10-241  10-242  10-243  10-244  10-245 	10-260  10-261  10-262  10-263  10-264  10-265 	10-350  10-351  10-352  10-353  10-354  10-355 	10-370  10-371  10-372  10-373  10-374  10-375 
465 X 280 X 135 mm 	10-246  10-247  10-248  10-249  10-250  10-251 	10-266  10-267  10-268  10-269  10-270  10-271 	10-356  10-357  10-358  10-359  10-360  10-361 	10-376  10-377  10-378  10-379  10-380  10-381 
465 X 280 X 150 mm 	10-252  10-253  10-254  10-255  10-256  10-257 	10-272  10-273  10-274  10-275  10-276  10-277 	10-362  10-363  10-364  10-365  10-366  10-367 	10-382  10-383  10-384  10-385  10-386  10-387 
Ps: All container measurements cover the outer size of the lids.				
				

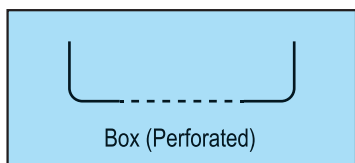

CODE	OUTSIDE	INSIDE	WEIGHT
10-300	475 X 275 X 95 mm	422 X 260 X 90 mm	1.65 kg
10-301	475 X 275 X 116 mm	422 X 260 X 96 mm	1.80 kg
10-302	475 X 275 X 137 mm	422 X 260 X 117 mm	1.95 kg

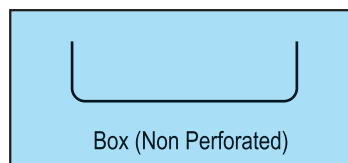

CODE	OUTSIDE	INSIDE	WEIGHT
10-305	475 X 275 X 95 mm	422 X 260 X 90 mm	1.50 kg
10-306	475 X 275 X 116 mm	422 X 260 X 111 mm	1.65 kg
10-307	475 X 275 X 137 mm	422 X 260 X 132 mm	1.75 kg



Perforated Lids



Safety Lids

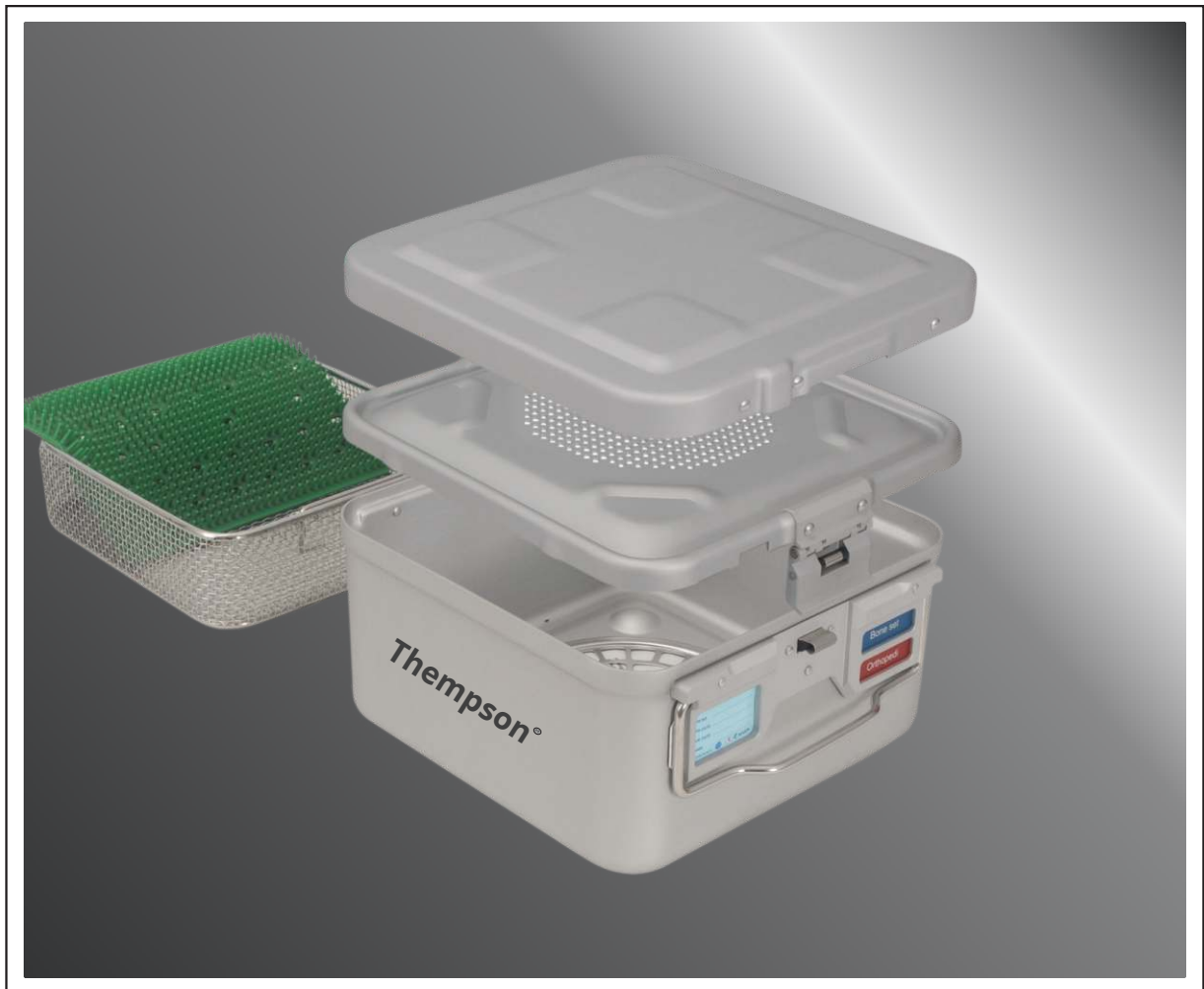


Removal Lids (Non Perforated)










Sterile Container System



1/2 Half Size

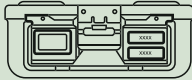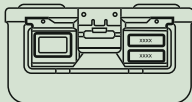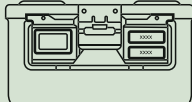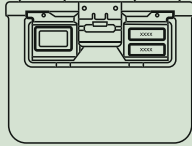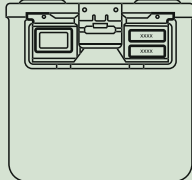Classic Model

Specifications Sizes	Lid (Perforated)	Lid (Perforated)	Safety Lid	Safety Lid
	Box (Non-Perforated)	Box (Perforated)	Box (Non-Perforated)	Box (Perforated)
285 X 280 X 100 mm 	10-500  10-501  10-502  10-503  10-504  10-505 	10-530  10-531  10-532  10-533  10-534  10-535 	10-640  10-641  10-642  10-643  10-644  10-645 	10-670  10-671  10-672  10-673  10-674  10-675 
285 X 280 X 135 mm 	10-506  10-507  10-508  10-509  10-510  10-511 	10-536  10-537  10-538  10-539  10-540  10-541 	10-646  10-647  10-648  10-649  10-650  10-651 	10-676  10-677  10-678  10-679  10-680  10-681 
285 X 280 X 150 mm 	10-512  10-513  10-514  10-515  10-516  10-517 	10-542  10-543  10-544  10-545  10-546  10-547 	10-652  10-653  10-654  10-655  10-656  10-657 	10-682  10-683  10-684  10-685  10-686  10-687 
285 X 280 X 200 mm 	10-518  10-519  10-520  10-521  10-522  10-523 	10-548  10-549  10-550  10-551  10-552  10-553 	10-658  10-659  10-660  10-661  10-662  10-663 	10-688  10-689  10-690  10-691  10-692  10-693 
285 X 280 X 260 mm 	10-524  10-525  10-526  10-527  10-528  10-529 	10-554  10-555  10-556  10-557  10-558  10-559 	10-664  10-665  10-666  10-667  10-668  10-669 	10-694  10-695  10-696  10-697  10-698  10-699 
Ps: All container measurements cover the outer size of the lids.				
				